SDR® Plus

Bulk Fill Flowable

Shade	Curing Time	
	2 mm	4 mm
Universal	20"	20"
A1 (pink)	20"	40"
A2 (black)	20"	40"
A3 (red)	20"	40"


Caution, consult
instructions for use

THE DENTAL
SOLUTIONS
COMPANY™


Dentsply
Sirona

Step-by-step instructions

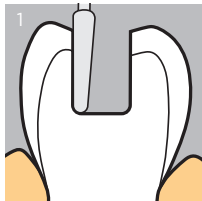
Cavity Preparation

Keep the cavity preparation to a minimum.

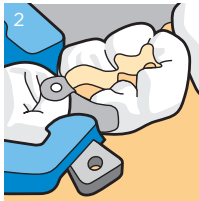
Matrix Placement

For Class II restorations place a matrix and wedge. Use of pre-contoured matrix bands is recommended.

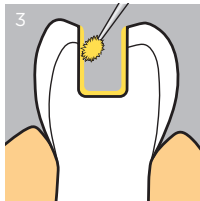Adhesive Application


Refer to adhesive manufacturer's directions for properly applying adhesive. Light cure taking depth of cavity into account.

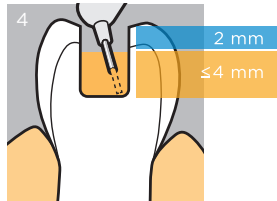SDR® Plus Application

Begin at the deepest portion of the cavity. Gradually withdraw while filling.


Step-by-step instructions

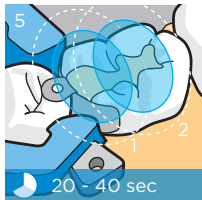Light Curing


Light cure with a minimum output of 550 mW/cm^2 .

U shade $\geq 20 \text{ sec}$.
A shades $\geq 40 \text{ sec}$.

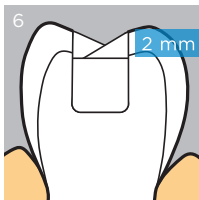Placement of Universal Composite

Restore occlusal 2 mm with oblique layers to facilitate occlusal contouring. Light cure according to the manufacturer's instructions for use.

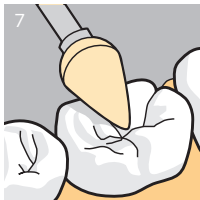Finishing & Polishing

Contour the restoration using finishing burs or diamonds. Finish the restoration by applying medium pressure and slowly back off to a light pressure.

For a glass-like luster, polish the restoration.

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SDR® Plus

Bulk Fill Flowable

CAUTION: For dental use only.
USA: Rx only

Content	Page
1 Product description	2
2 Safety notes	3
3 Step-by-step instructions	5
4 Hygiene	8
5 Lot number and expiration date	9

1 Product description

SDR® Plus Bulk Fill Flowable is a one-component, fluoride-containing, visible light cured, radio-paque resin composite restorative material. It is designed to be used as a base in Class I and II restorations. It is also suitable as a stand-alone restorative material in non-occlusal-contact applications.

SDR® Plus material has handling characteristics typical of a “flowable” composite, but can be placed in 4 mm increments with minimal polymerization stress.

SDR® Plus material has a self-leveling feature that allows intimate adaptation to the prepared cavity walls. When used as a base/liner, it is designed to be overlaid with a methacrylate based universal/posterior composite for replacing missing occlusal/facial enamel.

1.1 Indications

SDR® Plus material direct restorative:

- Base in cavity Class I & II direct restorations
- Liner under direct restorative materials – Class II box liner
- Pit & Fissure Sealant
- Conservative Class I restorations
- Core build-up
- Class III and V restorations

1.2 Contraindications

SDR® Plus material is contraindicated for use with patients who have a known hypersensitivity to methacrylate resins.

1.3 Delivery forms¹

SDR® Plus material is available in:

- Predosed Compula® Tips for direct intraoral application
- Small syringes for direct intraoral application
- One universal (approximates VITA®² B1 shade tab) and 3 VITA® (A1, A2, A3) shades

¹ Some delivery forms may not be available in all countries.

² Not a registered trademark of Dentsply Sirona, Inc.

1.4 Composition

- Barium-alumino-fluoro-borosilicate glass
- Strontium alumino-fluoro-silicate glass
- Ytterbium trifluoride glass
- Silicon dioxide
- Modified urethane dimethacrylate resin
- Polymerizable dimethacrylate resin
- Polymerizable trimethacrylate resin
- Triethyleneglycol dimethacrylate
- Camphorquinone (CQ) Photoinitiator
- Photoaccelerator
- Butylated hydroxyl toluene (BHT)
- UV Stabilizer
- Titanium dioxide
- Inorganic iron oxide
- Fluorescing agent
- Particles of inorganic filler range from 20 nm to 10 µm, total filler 47.3% by volume.

1.5 Compatible adhesives

SDR® Plus material is used following application of a suitable dentin/enamel adhesive and is chemically compatible with conventional methacrylate-based dentin/enamel adhesives including Dentsply Sirona adhesives designed for use with visible light cured composite restoratives (see complete Instructions for Use of selected adhesive).

1.6 Compatible enamel replacement restoratives

SDR® Plus material is used in conjunction with a suitable universal/posterior restorative material³ as an occlusal/facial enamel replacement (see Step-by-step instructions) and is chemically compatible with conventional methacrylate-based composite restorative materials including Dentsply Sirona visible light cured universal/posterior composite restorative materials designed for posterior Class I & II occlusal restorations (see complete Instructions for Use of selected occlusal restorative material).

2 Safety notes

Be aware of the following general safety notes and the special safety notes in other sections of these Instructions for Use.

Safety alert symbol.

- This is the safety alert symbol. It is used to alert you to potential personal injury hazards.
- Obey all safety messages that follow this symbol to avoid possible injury.

2.1 Warnings

SDR® Plus material contains polymerizable methacrylate monomers which may be irritating to skin, eyes and oral mucosa and may cause allergic contact dermatitis in susceptible persons.

- **Avoid eye contact** to prevent irritation and possible corneal damage. In case of contact with eyes, rinse immediately with plenty of water and seek medical attention.
- **Avoid skin contact** to prevent irritation and possible allergic response. In case of contact, red-dish rashes may be seen on the skin. If contact with skin occurs immediately remove material with cotton and alcohol and wash well thoroughly with soap and water. In case of skin sensitization or rash, discontinue use and seek medical attention.

³ Enamel replacement material overlay is not required for Pit & Fissure Sealant, Conservative Class I, Class III, Class V or core build-up applications.

- **Avoid contact with oral soft tissues/mucosa** to prevent inflammation. If accidental contact occurs, immediately remove material from the tissues. Flush mucosa with plenty of water after the restoration is completed and expectorate/evacuate the water. If inflammation of mucosa persists, seek medical attention.

2.2 Precautions

This product is intended to be used as specifically outlined in the Instructions for Use. Any use of this product inconsistent with the Instructions for Use is at the discretion and sole responsibility of the practitioner.

- Contact with saliva and blood during composite placement may cause failure of the restoration. Use of rubber dam or adequate isolation is recommended.
- Wear suitable protective eyewear, mask, clothing and gloves. Protective eyewear is recommended for patients.
- Devices marked "single use" on the labeling are intended for single use only. Discard after use. Do not reuse in other patients in order to prevent cross-contamination.
- The syringe cannot be reprocessed. To prevent the syringe from exposure to spatter or spray of body fluids or contaminated hands it is mandatory that the syringe is handled with clean/disinfected gloves. As additional precautionary measure, the syringe may be protected from gross debris but not from all contamination by applying a protective barrier.
- SDR® Plus material should extrude easily. **DO NOT USE EXCESSIVE FORCE.** Excessive pressure may result in unanticipated extrusion of the material or cause the Compula® Tip to eject from the Compules® Tips Gun.
- Use of Compula® Tips with the Dentsply Sirona Compules® Tips Gun is recommended.
- Tightly close syringes with original cap immediately after use. Do not allow dispensing tip to remain as a syringe cap.
- Only use the dispensing tips supplied by Dentsply Sirona for placement of the material.
- SDR® Plus material is available in semi-translucent shades. If placed to the cavosurface margin, e.g. occlusal, facial or lingual surface or wide proximal box extension cavities, the tooth/restorative demarcation may be visible. Heavily stained dentin may be visible through the cured SDR® Plus material. Use of an appropriately shaded and/or opaque restorative material of adequate thickness in visible areas/surfaces is recommended where esthetics is paramount.
- When placing SDR® Plus material as a base or liner, allow at least 2 mm for selected occlusal restorative material.
- Insufficient data exist to support the use as a Class I restoration with isthmus width beyond central groove areas, or when replacing a functional cusp. Excessive wear or restoration failure may result.
- Interactions:
 - Eugenol- and hydrogen peroxide-containing materials should not be used in conjunction with this product since they may interfere with hardening and cause softening of the polymeric components of the material.
 - SDR® Plus material is light cured material. Therefore, it should be protected from ambient light. Proceed immediately once the material has been placed.
 - If mineral-impregnated (e.g. ferric compounds) retraction cords and/or hemostatic solutions are used in conjunction with adhesive procedures, marginal seal may be adversely affected, allowing microleakage, subsurface staining and/or restoration failure. If gingival retraction is necessary, use of plain, non-impregnated cord is recommended.

2.3 Adverse reactions

Product may irritate the eyes and skin.

- Eye contact: Irritation and possible corneal damage.
- Skin contact: Irritation or possible allergic response. Reddish rashes may be seen on the skin.
- Mucous membranes: Inflammation (see Warnings).

2.4 Storage conditions

Inadequate storage conditions will shorten the shelf life and may lead to malfunction of the product.

- Keep out of direct sunlight and store in a well ventilated place at temperatures between 2-24 °C/35-75 °F.
- Allow material to reach room temperature prior to use.
- Protect from moisture.
- Do not freeze.
- Do not use after expiration date.

3 Step-by-step instructions

3.1 Cavity preparation

1. Prepare the cavity so that no residual amalgam or restorative material is left.
2. Rinse surface with water spray and carefully dry it with air spray. Do not desiccate the tooth structure.
3. Use a dental dam or cotton rolls to isolate the cavity from contamination.

3.2 Placement of matrix

For optimal proximal contacts proceed as follows:

1. Place a matrix (e.g. AutoMatrix® matrix system or Palodent® Plus/Palodent® V3 Sectional matrix system) and wedge. Burnishing of the matrix band will improve contact and contour. Pre-wedging/ ring placement is advocated.

3.3 Pulp protection, tooth conditioning/dentin pretreatment, adhesive application

Refer to adhesive manufacturer's directions for pulp protection, tooth conditioning and/or adhesive application. Once the surfaces have been properly treated, they must be kept uncontaminated. Proceed immediately to placement of SDR® Plus material.

3.4 Dispensing SDR® Plus material

Danger of injury due to excessive force.

1. Apply slow and steady pressure on the syringe plunger or applicator gun.
2. Do not use excessive force – Syringe or Compula® Tip rupture or ejection from applicator gun may result.

Compula® Tips

The pre-dosed Compula® Tip provides the combination of a unit dose Compules® Tip with a metal cannula Applicator Tip.

1. Load Compules® Tips Gun with pre-dosed Compula® Tip. Insert Compula® Tip into the notched opening of the Compules® Tips Gun barrel. Be certain that the collar on the Compula® Tip is inserted first.
2. Remove the colored cap from the Compula® Tip. The Compula® Tip may be rotated 360° to gain the proper angle of entrance into the cavity.
3. Dispense the material into the cavity preparation using a slow, steady pressure. **DO NOT USE EXCESSIVE FORCE.**
4. To remove the used Compula® Tip, be sure that the Compules® Tips Gun plunger is pulled back completely by allowing the handle to open to its widest position. Apply a downward motion to the front end of the Compula® Tip and remove.

Syringes

1. Remove cap from the end of the syringe. To assure free flow of material from syringe, express a small amount onto pad, away from the patient field.
2. Attach disposable, black dispensing tip to end of the syringe. Turn tip clockwise $\frac{1}{4}$ to $\frac{1}{2}$ turn to assure that it is fully seated. Tug on tip to be sure that it is locked into the collar of the syringe.
3. Material should flow freely with gentle pressure. DO NOT USE EXCESSIVE FORCE. If more than gentle pressure is required, remove from patient field and check for obstruction.
4. DISCARD AND PROPERLY DISPOSE DISPENSING TIP IMMEDIATELY AFTER USE. REPLACE ORIGINAL CAP. Do not store syringe with black dispensing tip in place. STORE ONLY WITH ORIGINAL CAP.

Note: It is recommended to pull back slightly on the syringe plunger after use to prevent excessive flow of material.

3.5 Placement of SDR® Plus material

Note: As SDR® Plus material is designed to conform to cavity shape, proper contact areas must be established with the matrix system selected. Material will not distend matrix band. If needed, manually distend matrix band with a suitable hand instrument such as a condenser or ball burrisher prior to and/or during light curing.

1. Dispense SDR® Plus material directly into preparation site/tooth surface using slow, steady pressure. Begin dispensing at the deepest portion of the cavity, keeping tip close to cavity floor. Gradually withdraw tip as cavity is filled. Avoid lifting the tip out of dispensed material while dispensing to minimize air entrapment. At the completion of dispensing, wipe tip against cavity wall while withdrawing from the operative field.
2. Within a few seconds, dispensed SDR® Plus material will self-level, eliminating the need for further manipulation with hand instruments. In case of overfill, use a flocked applicator tip wetted with residual adhesive to remove excess. Any visible air bubbles should be pierced with a clean, sharp explorer prior to curing.
3. When used as a Bulk Fill Base material, most cavities may be filled in one bulk increment (up to 4 mm) as needed to fill the cavity 2 mm short of the occlusal cavosurface⁴. In deeper preparations, place material in 4 mm increments, thoroughly light curing each increment.
4. Alternatively, SDR® Plus material may be placed in a thin layer as a traditional "flowable" liner on exposed dentin.

3.6 Curing

SDR® Plus material is designed to be cured in increments up to a 4 mm depth/thickness.

1. Light cure each area of the restoration surface with a suitable visible light curing unit designed to cure materials containing camphorquinone (CQ) initiator, i.e. spectral output containing 470 nm. Minimum light output must be at least 550 mW/cm².

Pulp Injury due to excessive heat – curing output > 2,000 mW/cm².

1. Check Curing Light manufacturer's literature for stated output.
2. Consult Curing Light Manufacturer's Directions for compatibility curing recommendations.
3. Do not use table below for recommended curing times with output > 2,000 mW/cm².

⁴ When used as a Pit & Fissure Sealant, conservative Class I, Class III or V restorative or core build-up material, material may be placed to full contour including up to the cavosurface margin.

Curing Recommendations ⁵			
Shade	Light Output	Cure Time 2 mm	Cure Time 4 mm
Universal	Halogen and LED Lights 550-1,000 mW/cm ²	20 sec.	20 sec.
	High Power LED Lights 1,000-2,000 mW/cm ²	10 sec.	10 sec.
A1, A2, A3	Halogen and LED Lights 550-1,000 mW/cm ²	20 sec.	40 sec.
	High Power LED Lights 1,000-2,000 mW/cm ²	10 sec.	25 sec.

Inadequate polymerization due to insufficient curing.

1. Check compatibility of curing light.
2. Check curing cycle.
3. Check curing output before each procedure.

3.7 Completion

1. When used as a base/liner, proceed immediately to complete the restoration using a methacrylate based universal/posterior restorative material following manufacturer's Instructions for Use.

Technique Tip: In most applications, a body shade of occlusal restorative material provides the best esthetic result. In case of heavily stained dentin, use of an opaque shade may be required. Do not disturb or allow contamination of cured surface or exposed adhesive. If contamination occurs, follow adhesive manufacturer's Instructions for Use for re-application of adhesive. If surfaces remain undisturbed, additional application of adhesive between increments is neither necessary nor recommended.

3.8 Finishing and polishing

1. If contouring, finishing, and/or polishing of cured SDR[®] Plus material is necessary, use traditional rotary devices following manufacturer's Instructions for Use.
2. Removal of resin flash is best accomplished with the Enhance[®] Finishing System. The Enhance[®] System will remove flash and finish restoration margins.
3. Polish final restoration using Enhance[®] PoGo[®] Polishing system or Prisma[®] Gloss[™] Polishing Paste. (See complete Instructions for Use supplied with polishing product chosen).

All shades of SDR[®] Plus material are radio-opaque, with a 1 mm radio-opacity equivalent to 2.6 mm radio-opacity of aluminum. Aluminum has a radio-opacity equivalent to that of dentin. Thus 1 mm of material having a radio-opacity equivalent to 1 mm of aluminum has a radio-opacity equivalent to that of dentin.

⁵ Refer to curing light manufacturer's recommendations for compatibility and curing recommendations.

4 Hygiene

Cross-contamination.

Infection.

1. Do not reuse single use products. Dispose of in accordance with local regulations.
2. Reprocess reusable products according to the instructions.

4.1 Syringes – cross-contamination

Cross-contamination.

Infection.

1. The syringes cannot be reprocessed.
2. To prevent the syringes from exposure to spatter or spray of body fluids or contaminated hands it is mandatory that the syringes are handled with clean/disinfected gloves. Do not reuse syringes if contaminated.
3. Dispose of contaminated syringes in accordance with local regulations.

For syringe, remove used dispensing tip and discard appropriately. Replace original syringe cap prior to storing. To prevent syringes from exposure to spatter or spray of body fluids or contaminated hands, or oral tissues, use of a protective barrier is recommended. The use of protective barriers is an additional precautionary measure against gross debris but not against all contamination.

Incidental contact of the syringe with water, soap or a water-based hospital-level disinfection solution will not damage syringe body. Do not allow any solution contact with contained material. Discard composite material that has been in contact with any fluid or non-sterile instrument.

Repeated liquid contact may damage label. Dry the syringe with a lint-free single use cloth.*

*** Note: Destruction of the label by vigorous wiping.**

Wipe syringe gently.

4.2 Compules® Tips Gun

For reprocessing instructions please refer to the Instructions for Use of the Compules® Tips Gun, which is available on our webpage at www.dentsplysirona.com/IFU. If requested, we will send you a free printed copy of the Instructions for Use in the language you require within 7 days. In the USA, call 1-800-532-2855. Outside North America, use the order form provided on www.dentsplysirona.com/IFU for this purpose.

4.3 Compula® Tips – cross-contamination

Cross-contamination.

Infection.

1. Do not reuse Compula® Tips. Dispose Compula® Tips in accordance with local regulations.

5 Lot number () and expiration date ()

1. Do not use after expiration date.
ISO standard is used: "YYYY-MM" or "YYYY-MM-DD".
2. The following numbers should be quoted in all correspondence:
 - Reorder number
 - Lot number
 - Expiration date

© 2018 Dentsply Sirona, Inc. (2018-03-30). All Rights Reserved.

[These Instructions for Use are based on Master Version 03]

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© 2018 Dentsply Sirona, Inc. (2018-03-30). Todos os direitos reservados.
Dentsply DeTrey GmbH