

## Lingual Attachments



*Our Strength is our Service*

# Lingual Attachments

## Preformed Bondable Lingual Retainers

- Central to central
- Cuspid to cuspid
- Bicuspid to bicuspid


Universal shape - can be used to retain both upper and lower cuspid - cuspid and bicuspid - bicuspid. Also, one universal size helps keep inventory to a minimum. Sizes are permanently laser etched for quick and easy identification. Sold in a kit, or individually.

### Lingual Retainer Starter Kit

Kit includes 8 popular sized retainers (3-I-3) **O/No 2165-200**  
 1 central-central retainer (1-I-1)

### Individual Lingual Retainers

Size	O/No	Size	O/No
10mm	<b>2165-010</b>	33mm	<b>2165-330</b>
22mm	<b>2165-220</b>	34mm	<b>2165-340</b>
23mm	<b>2165-230</b>	35mm	<b>2165-350</b>
24mm	<b>2165-240</b>	36mm	<b>2165-360</b>
25mm	<b>2165-250</b>	37mm	<b>2165-370</b>
26mm	<b>2165-260</b>	39mm	<b>2165-390</b>
27mm	<b>2165-270</b>	40mm	<b>2165-400</b>
28mm	<b>2165-280</b>	41mm	<b>2165-410</b>
29mm	<b>2165-290</b>	42mm	<b>2165-420</b>
30mm	<b>2165-300</b>	43mm	<b>2165-430</b>
31mm	<b>2165-310</b>	44mm	<b>2165-440</b>
32mm	<b>2165-320</b>		

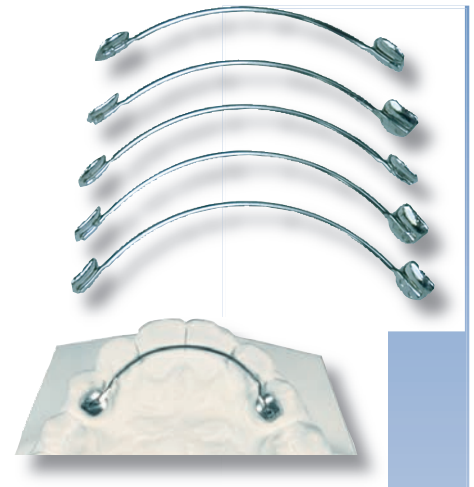

## .0195 Lingual Retainer Wire on Spool

Fabricate retainers by bonding a length of heavy multistrand wire to the lingual surfaces of anterior teeth, to act as a 'splint' after removal of orthodontic bands/brackets

- Six strand coaxial wire
- Economical 30 foot spool
- No wastage

We recommend using the Wire Bonder with this type of retainer - see following page

### .0195 Lingual Retainer Wire on Spool


30ft/9.144m per reel **O/No 5300-395**



## Weldable Stripform Lingual Cleats

### Weldable Stripform Lingual Cleats

100 per pack **O/No 5028-317**

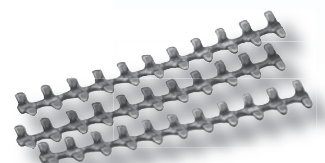

## The 'Button Bonder'



The Button Bonder allows you to make tooth coloured lingual buttons using *your own light cured composite*. Simply fill the hollow (button mould) tip with light cure composite, hold against the tooth and light cure. Remove the tip and you now have a functional lingual button at a fraction of the cost of a metal button or cleat (Smoother and more rounded too!).

Each Button Bonder tip can be cold sterilised and re-used.

### Button Bonder

10 per pack	O/No <b>9000-609</b>
-------------	-------------------------

## The Wire Bonder

It's Quick... It's Easy... It Works!

Helps you to produce neat lingual retainers without 'the mess'! The Wire Bonder is a tiny translucent 'cookie cutter' with a built-in wire slot. Fill the tip with light cure resin, place on the tooth over the wire, put the wire into the slot, and hold in desired location on the tooth and light cure. You now have a bonded wire with a small, precise, smooth mound of resin holding the wire to the tooth. Lingual retainers can be made in a few minutes.

- Quick, easy and precise - no rough edges
- Forms a small precise mound of light-cure resin bonding the wire to the tooth
- Smooth, neat composite attracts less plaque and is easier to clean
- Uniform thickness covering the wire
- Re-use over and over - Cold Sterilise
- Can be used on the labial or the lingual surfaces
- Alignment tabs make it easy to position the tip correctly on the teeth

### Wire Bonder

10 tips and handle per pack	O/No <b>9000-608</b>
-----------------------------	-------------------------

We recommend the use of 'Solar' Light Cure Adhesive (see page XXXX) for the best results with the Wire Bonder one minute bonded retainers. For the full 'Bonder Line' see page XXXX.

## Weldable High Profile Buttons

### Weldable High Profile Button

Curved	O/No <b>5028-370</b>
Flat	<b>5028-360</b>

## Weldable Lingual Sheaths

### Weldable Lingual Sheaths

10 per pack	O/No
For use with removable Quad Helix	<b>5000-314</b>
Curved length 3.5mm	<b>5000-314</b>
For use with palatal arches	<b>5000-320</b>



# Lingual Attachments

## Direct Bond Lingual Cleats

On mini mesh base. Universal contour.

### Lingual Attachments

10 per pack

O/No  
**5028-170**



## Direct Bond High Profile Button

On mini mesh base.

### Direct Bond High Profile Button

10 per pack

Flat

Curved

O/No  
**5028-160**  
**5028-165**



## Direct Bond Eyelets

On mini mesh base. Universal contour.

### Direct Bond Eyelets

10 per pack

O/No  
**5028-150**



## 'Technique' Cuspid Traction System

14 K Gold Chain and Pad NO ETCHING REQUIRED

Gold chain pre-attached to a bonding base.

'Technique' provides an easy way to attach to a pre-emerged tooth, so an orthodontic device can be used to draw the tooth into position.

Individual bonding pads for cuspids, bicuspids, molars and incisors.

'Technique' requires no etching. Just expose the tooth surface, dry and bond base in place using the glass ionomer adhesive.

### Technique Cuspid Traction System

Gold Chain and Pad 10 per pack	O/No
Cuspid	<b>5028-005</b>
Bicuspid	<b>5028-007</b>
Incisor	<b>5028-008</b>
Individual Gold Chain and Pad	
Cuspid	<b>5028-005I</b>
Bicuspid	<b>5028-007I</b>
Incisor	<b>5028-008I</b>
Technique cement powder 15gm	<b>5028-003</b>
Technique cement liquid 15cc	<b>5028-004</b>

