Zhermack
Dental

Zeta 7 Spray

EN - Instructions for use 2

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SV - Bruksanvisning 30

DA - Brugsanvisning 34

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LT - Naudojimo instrukcija 46

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ET - Kasutusjuhised 54

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SK - Návod na použitie 62

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HR - Upute za uporabu 74

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Zhermack
Dental

Zeta 7 Spray

Disinfectant and cleaner for dental impressions

MD


CE 0123

dentsplysirona.com/zhermack

Zhermack S.p.A. | Via Bovazecchino, 100 | 45021 Badia Polesine (RO) Italy - T +39 0425 597611 | F +39 0425 53596 | info@zhermack.com | www.zhermack.com

2100-EIFU-017-001 Zeta7_Spray

Last update: 2026-05

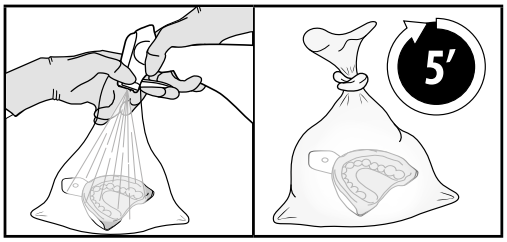1 2

EN – Zeta 7 Spray
Instructions for use

1. INTENDED USE

Disinfectant and cleaner for dental impressions.

2. PRODUCT DESCRIPTION

Zeta 7 Spray is a ready-to-use alcohol-based disinfectant for cleaning and disinfecting dental impressions.

3. PACKAGING

750 ml bottle.

4. COMPOSITION

100 g of Zeta 7 Spray contain 83.5 g of ethanol, 10 g of 2-propanol, non-ionic surfactants, additives, auxiliaries and water q.s. for 100 g.

5. INDICATIONS FOR USE

The device is intended for use in the dental sector, in professional dental surgeries (public or private health facilities with the necessary legal authorisations), by professionals trained and qualified in cleaning and disinfecting dental impressions.

6. CONTRAINDICATIONS

Use of the product is contraindicated for operators with known hypersensitivity to limonene and citral.

7. SIDE EFFECTS

Allergic reactions may occur following contact with Zeta 7 Spray. Direct contact of the product with eyes causes serious eye irritation; contact with skin can cause skin irritation. Inhaling its vapours can cause respiratory tract irritation.

8. CLINICAL BENEFITS

Zeta 7 Spray helps reduce the microbial load of dental impressions after taking the impression, decreasing the risk of infection transmission, if used strictly in accordance with the instructions for use and the safety data sheet.

9. WARNINGS

- This product contains limonene and citral which may cause an allergic skin reaction. Use by operators sensitive to limonene and citral is contraindicated. If irritation, redness or other signs of hypersensitivity occur, stop using the product and seek medical advice.
- Harmful if ingested. Do not ingest the product. If ingested, DO NOT induce vomiting and do not give anything to eat or drink. In case of ingestion and/or contact, consult a poison centre or seek medical advice immediately. Avoid direct contact with the product.
- Causes serious eye irritation. Avoid direct contact with eyes. Protect eyes with airtight goggles. In case of contact with eyes, rinse immediately with plenty of water for several minutes and seek medical advice if irritation persists.
- May irritate the respiratory tract. Do not inhale vapours. Wear a face mask (see section 8 of the product safety data sheet), adequately ventilate places in which the product is handled and ensure that the bottle is stored tightly closed.
- May cause skin irritation. Avoid direct contact with the skin and face. Protect hands with gloves, wear a face mask, protective clothing and safety footwear for professional use (see section 8 of the product safety data sheet). In the event of accidental contact with the skin, wash thoroughly with soap and water and seek medical advice if necessary.
- Flammable liquid and vapour. Keep away from heat, sparks, naked flames, hot surfaces and other sources of ignition. Ensure the bottle is stored tightly closed in a well-ventilated place.
- Only use the product in the operating area when there are no patients and/or individuals not wearing appropriate PPE present.
- To reduce the risk of cross-contamination, always use new, uncontaminated gloves when handling the products (bottle and spray dispenser).

10. PRECAUTIONS

- Keep the product's instructions for use for the entire duration of its use.
- Zeta 7 Spray is compatible with the following impression materials: addition and condensation silicone, alginate, polyether, and polysulphide. For further limitations and indications on use regarding the disinfection of these materials, read the instructions for use of the devices to be disinfected.
- Zeta 7 Spray is compatible with both metal (stainless steel, aluminium) and plastic (polycarbonate, polystyrene) impression trays.
- This disinfectant was not designed for reprocessing impression trays. To reprocess impression trays, follow the manufacturer's instructions for use.

• Avoid contact of Zeta 7 Spray liquid with other disinfectants or cleaners.

• Zeta 7 Spray must be used at temperatures above 20 °C. If frozen, it must be thawed and maintained at room temperature for at least 24 hours before use.

• Always comply with the contact times indicated in section 11. Exposure for longer times may result in deterioration of the treated surfaces.

• Do not use the product and dispose of it in accordance with the instructions and the product safety data sheet:

- if the batch, expiry and UDI information on the outer packaging is missing or illegible;
- if the primary packaging (bottle, spray dispenser) is damaged or does not work properly.

• Do not use after the use-by date shown on the package, indicated on the front of the label in year-month-day format.

• Once the contents of the bottle are finished, dispose of the packaging as indicated in section 14.

• Do not use the bottle for purposes other than those indicated in these instructions for use.

11. STEP-BY-STEP INSTRUCTIONS

11.1 PRELIMINARY OPERATIONS

1.Read the instructions for use and the product safety data sheet before use.

2.Adequately ventilate the rooms where the product is handled.

3.Protect your hands with a pair of gloves, protect your eyes and face with a pair of airtight goggles and wear a face mask, work clothes and safety footwear for professional use (see section 8 of the product safety data sheet).

11.2 MANUAL CLEANING OF IMPRESSIONS

1. Rinse the impression under tap water for **60 seconds** immediately after removal from the mouth. Clean the impression within 5 minutes of removal from the mouth. If necessary, increase the rinsing time until no dirt remains visible.

2. Remove any excess water from the impression (eliminate any visible water residue).

11.3 MANUAL DISINFECTION OF IMPRESSIONS

1. After cleaning (see section 11.2), place the dental impression in a transparent bag.

2. Insert the nozzle of the Zeta 7 Spray bottle into the transparent bag and spray **10 times** to cover the entire surface of the dental impression and the impression tray (fig.1).

3. Hold the bag closed to prevent aerosol escape.

4. Leave the dental impression in the transparent bag for **5 minutes** of contact with Zeta 7 Spray, keeping the bag closed (fig.2).

5. Remove the impression from the bag and dispose of the bag.

6. If required by the manufacturer's instructions for the impression material being disinfected, rinse and/or dry.

12. INFORMATION ON PROPER MANAGEMENT OF PACKAGING

The product packaging can be safely reused without undergoing cleaning and disinfection processes, provided that the following storage and handling recommendations are followed:

• Use the product only when the patient is outside the operating area;

• Always wear new, uncontaminated gloves when handling packaging. If the gloves get contaminated, replace them immediately before handling the product;

• Avoid placing packaging in contact with contaminated surfaces and/or devices;

• Store the tightly closed package outside the operating area in a clean, closed and dry environment protected from splashes, sprays or aerosols of body fluids.

In case of contamination, dispose of packaging and its content as indicated in section 14 dedicated to disposal.

13. SPECTRUM OF ACTION

Table 1

Effectiveness (micro-organisms)	**Reference standard** *Test carried out in clean conditions*
Bactericidal (S. aureus, P. aeruginosa, E. hirae)	EN 13727, EN 17387
Yeasticidal (C. albicans)	EN 13624, EN 17387
Mycobactericidal including tuberculocidal (M. terrae, M. avium)	EN 14348
Virucidal (poliovirus, adenovirus, rotavirus, norovirus, including HIV, HBV, HCV)	EN 14476, EN 16777

Activity in suspension is tested in compliance with the standards EN 13727, EN 13624, EN 14348, EN 14476. Surface activities are tested in compliance with the standards EN 17387 and EN 16777.

14. STORAGE, STABILITY AND DISPOSAL

Store the product tightly closed outside the operating area in a clean, closed, cool and dry environment protected from splashes, sprays or aerosols, at a maximum temperature of 27°C/80°F, away from heat sources, sparks, open flames and hot surfaces and other sources of ignition and direct sunlight. Do not use after the use-by date indicated on the packaging.

The following disposal instructions must be considered valid only if the product has been stored and used according to these instructions for use:

• Waste with only chemical residues (residues of disinfectant solution) are a physical hazard for people and the environment and must be handled as **hazardous waste** in accordance with local regulations.

• Waste with biological residues must be disposed of as **special waste** at risk of biological contamination.

15. IMPORTANT REMARKS

Information provided in any way, even during demos, does not invalidate the instructions for use. Operators are required to check that the product is suitable for the application envisaged. The manufacturer cannot be held liable for damage, including to third parties, deriving from failure to follow instructions or from unsuitability for an application. The manufacturer's liability is, in any case, limited to the value of the products supplied. Report any serious incident involving the medical device to the manufacturer and to the relevant authorities.

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